

Importance of Microcystins/Nodularins Determination

Most of the world's population relies on surface freshwaters as its primary source for drinking water. The drinking water industry is constantly challenged with surface water contaminants that must be removed to protect human health. Toxic cyanobacterial blooms are an emerging issue worldwide due to increased source water nutrient pollution caused by eutrophication. Microcystins and Nodularins are cyclic toxin peptides. Microcystins (of which there are many structural variants or congeners) have been found in fresh water throughout the world. To date, approximately 140 variants of Microcystin have been isolated. The most common variant is Microcystin-LR. Other common Microcystin variants include YR, RR, and LA. These toxins are produced by many types of cyanobacteria (blue-green algae), including *Microcystis*, *Anabaena*, *Oscillatoria*, *Nostoc*, *Anabaenopsis*, and terrestrial *Hapalosiphon*. Nodularins are produced by the genus *Nodularia* and are found in marine and brackish water.

Acute poisoning of humans and animals constitutes the most obvious problem from toxic cyanobacterial blooms, and in several cases has led to death. Human and animal exposure to these toxins occurs most frequently through ingestion of water, through drinking or during recreational activities in which water is swallowed. These toxins mediate their toxicity by inhibiting liver function and are potent inhibitors of the serine/threonine protein phosphatases, and therefore may act as tumor promoters.



To protect consumers from adverse health effects caused by these toxins, the World Health Organization (WHO) has proposed a provisional upper limit for Microcystin-LR of 1.0 ppb ($\mu\text{g/L}$) in drinking water. The U.S. Environmental Protection Agency (EPA) has also established guidelines for Microcystins in drinking water at 0.3 ppb ($\mu\text{g/L}$) for children below school age and 1.6 ppb ($\mu\text{g/L}$) for all other age groups.

Performance Data

Test sensitivity: The detection limit for this assay (90% B/B₀ calculated from the average of 46 standard curves), based on MC-LR, is 0.016 ppb ($\mu\text{g/L}$) (see representative curve, below left).

Test reproducibility: The standard curve R² must be ≥ 0.98 . The absorbance Coefficient of variation (CVs) for standards should be $\leq 10\%$ and for the samples and Control should be $\leq 15\%$. The Control should be within its acceptable range and Standard 0 absorbance value should be between 0.8 - 3.000.

Selectivity*: The antibody exhibits very good cross-reactivity with all cyanobacterial cyclic peptide toxin congeners tested to date (see cross-reactivity illustration, below right).

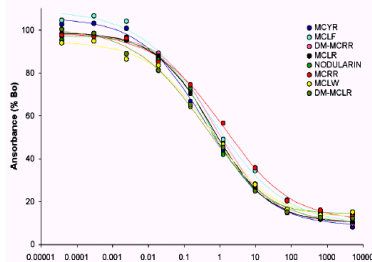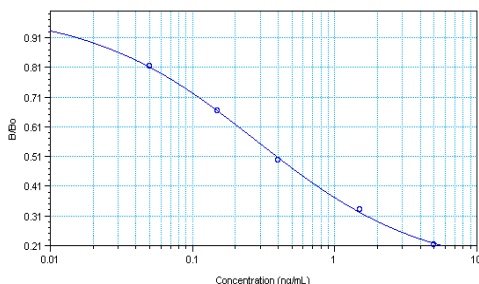

For demonstration purposes only. Not for use in sample interpretation.

References

- (1) W. J. Fischer, I. Garthwaite, C.O. Miles, K.M. Ross, J.B. Aggen, A.R. Chamberlin, N.A. Towers, and D.R. Dietrich, Congener-Independent Immunoassay for Microcystins and Nodularins. *Environ. Sci. Technol.* 35, 2001, 4849-4858.
- (2) Worldwide Patenting PCT WO 01/18059 A2.
- (3) U.S. Patent Number 6,967,240.
- (4) U.S. Patent 9,739,777.

†QuikLyse™ reagents may be used in a method of U.S. Patent 9,739,777

General Limited Warranty/Disclaimer: Gold Standard Diagnostics warrants the products manufactured by the Company, against defects and workmanship when used in accordance with the applicable instructions for a period not to extend beyond the product's printed expiration date. **Gold Standard Diagnostics makes no other warranty, expressed or implied. There is no warranty of merchantability or fitness for a particular purpose.**

For ordering or technical assistance contact:

Gold Standard Diagnostics

795 Horsham Road

Horsham, PA 19044

WEB: www.abraxiskits.com

Date this User Guide is effective: 12/11/2024

Tel: (215) 357-3911

Fax: (215) 357-5232

Ordering: info.abraxiskits@us.goldstandarddiagnostics.com

Technical Support: support.abraxiskits@us.goldstandarddiagnostics.com

Version: 04



ABRAXIS® Microcystins-ADDA SAES ELISA Microtiter Plate

Streptavidin Amplified Enhanced Sensitivity Enzyme-Linked
Immunosorbent Assay for the Congener-Independent* Determination of
Microcystins and Nodularins in Water Samples
Product No. 520011SAES

1. General Description

The ABRAXIS® Microcystins-ADDA SAES ELISA is an enhanced sensitivity immunoassay for the quantitative and sensitive congener-independent* detection of Microcystins and Nodularins in water samples. This test is suitable for the quantitative and/or qualitative detection of Microcystins and Nodularins in water samples. Please see the Gold Standard Diagnostics website for the appropriate technical bulletins for sample collection, handling, treatment of drinking (treated and untreated) and recreational water samples. If necessary, positive samples can be confirmed by HPLC, protein phosphatase assay, or other conventional methods.

2. Safety Instructions

The standard solutions in the test kit contain small amounts of Microcystins. The substrate solution contains tetramethylbenzidine (TMB) and the stop solution contains diluted sulfuric acid. Avoid contact of the TMB and stopping solution with skin and mucous membranes. If these reagents come in contact with skin, wash with water.

3. Storage and Stability

The ABRAXIS® Microcystins-ADDA SAES ELISA kit should be stored in the refrigerator (2-8°C). The solutions must be allowed to reach room temperature (20-25°C) before use. Reagents may be used until the last day of the month as indicated by the expiration date on the box. Consult state, local, and federal regulations for proper disposal of all reagents.

4. Microcystins-ADDA Antibody ELISA Test Principle

The test is an indirect competitive ELISA for the congener-independent detection of Microcystins and Nodularins. It is based on the recognition of Microcystins, Nodularins, and their congeners by specific biotinylated antibodies. Toxin, when present in a sample and a Microcystins-protein analogue immobilized on the plate compete for the binding sites of the biotinylated anti-Microcystins/Nodularins antibodies in solution. The plate is then washed and a streptavidin-HRP is added. After a second washing step and addition of the substrate solution, a color signal is generated. The intensity of the blue color is inversely proportional to the concentration of Microcystins present in the sample. The color reaction is stopped after a specified time and the color is evaluated using an ELISA reader. The concentrations of the samples are determined by interpolation using the standard curve constructed with each run.

5. Limitations of the ABRAXIS® Microcystins-ADDA SAES ELISA, Possible Test Interference

Numerous organic and inorganic compounds commonly found in water samples have been tested and found not to interfere with this test. However, due to the high variability of compounds that may be found in water samples, test interferences caused by matrix effects cannot be completely excluded.

Samples containing methanol must be diluted to a concentration of < 1% methanol to avoid matrix effects.

No matrix effects have been observed with seawater samples (salinity up to 38 parts per thousand).

No matrix effects have been observed with samples, which have been treated with sodium thiosulfate at concentrations ≤ 10 mg/mL or ascorbic acid at concentrations ≤ 1 mg/mL. *Note: Ascorbic acid used for dechlorination of treated drinking water samples may degrade Microcystins and may result in inaccurate (falsely low) results. See Sample Collection and Handling (section C) for information on the quenching of chlorinated water samples.*

Mistakes in handling the test can cause errors. Possible sources for such errors include: inadequate storage conditions of the test kit, incorrect pipetting sequence or inaccurate volumes of the reagents, too long or too short incubation times during the immune and/or substrate reaction, and extreme temperatures during the test performance (lower than 10°C or higher than 30°C).

As with any analytical technique (GC, HPLC, etc.), positive results requiring regulatory action should be confirmed by an alternative method.

6. Working Instructions

A. Materials Provided

1. Microtiter plate (12 X 8 strips) coated with an analog of Microcystins conjugated to a protein
2. Standards (6): 0, 0.05, 0.15, 0.4, 1.5, 5.0 ppb, 1.5 mL each
3. Control: 0.75 ± 0.185 ppb, 1.5 mL, prepared from a secondary source, for use as a Quality Control Standard (QCS)
4. Sample Diluent, 25 mL, for use as a Laboratory Reagent Blank (LRB) and for dilution of samples above the range of the standard curve
5. Microcystins-ADDA SAES Antibody Solution, 6 mL
6. Microcystins-ADDA SAES Conjugate Solution, 12 mL
7. ABRAXIS® Wash Buffer (5X) Concentrate, 100 mL, must be diluted prior to use, see Test Preparation (Section E)
8. Substrate (Color) Solution (TMB), 12 mL
9. Stop Solution, 6 mL

B. Additional Materials (not delivered with the test kit)

1. Micro-pipettes with disposable plastic tips (20-200 µL)
2. Multi-channel pipette (50-300 µL), stepper pipette (50-300 µL), or electronic repeating pipette with disposable plastic tips
3. Deionized or distilled water
4. Container with 500 mL capacity (for diluted 1X Wash Buffer, see Test Preparation, Section E)
5. Graduated cylinder
6. Paper towels or equivalent absorbent material
7. Timer
8. Tape or parafilm
9. Microtiter plate reader (wavelength 450 nm)
10. Microtiter plate washer (optional)

C. Sample Collection and Handling

Collect water samples in **glass or PETG** containers and test within 24 hours. Use of other types of plastic collection and/or storage containers may result in adsorptive loss of Microcystins, producing inaccurate (falsely low) results. Drinking water samples which contain chlorine or other oxidizers must be treated with sodium thiosulfate immediately after collection (refer to appropriate technical bulletin). If samples must be held for longer periods (up to 5 days), samples should be stored refrigerated. For storage periods greater than 5 days, samples should be stored frozen.

If total Microcystins concentration (free and cell bound) is required, an appropriate cell lysing procedure (freeze and thaw, ABRAXIS® QuikLyse™₊, etc.) must be performed prior to analysis. *Note: The use of sonication in cell lysing can negatively affect toxin concentrations, producing falsely low sample results. Please see the appropriate sample preparation technical bulletin for additional information on cell lysis.*

Samples may be filtered prior to analysis using glass fiber filters. The use of alternate filter types (non-glass fiber filters) may produce falsely low sample results, as Microcystins may bind to the filter material, removing it from the sample. Also, please note that some glass fiber filters are manufactured using a process, which may cause interference, which would cause inaccurate (falsely high) results. To avoid this potential bias in sample results, a total volume of at least 10 mL should be passed through the glass fiber filter, with the first 5 mL of filtered sample being discarded and the second 5 mL collected for testing (please see the Microcystins sample filtration technical bulletin for additional information on sample filtration). If determining total Microcystins concentration, samples should be lysed prior to filtration to prevent the removal of cell-bound Microcystins, which would cause inaccurate (falsely low) results.

D. Notes and Precautions

Micro-pipetting equipment and pipette tips for pipetting the standards and the samples are necessary.

The use of a multi-channel pipette, stepping pipette, or electronic repeating pipette is recommended for the addition of the antibody, enzyme conjugate, substrate, and stop solutions in order to equalize the incubation periods on the entire microtiter plate.

To avoid drift and obtain accurate results, the addition of the antibody, conjugate, color, and stop solutions should be performed in less than 2 minutes for each reagent. If additions to the entire microtiter plate cannot be completed in less than 2 minutes, run size should be decreased to the number of rows which can be pipetted in less than 2 minutes.

Please use only the reagents and standards from one kit lot in one test, as they have been adjusted in combination.

E. Test Preparation

1. Allow the reagents and samples to reach ambient temperature before use.
2. Remove the number of microtiter plate strips required from the resealable pouch. The remaining strips are stored in the pouch with the desiccant (tightly sealed).
3. The standards, control, sample diluent (LRB), antibody, enzyme conjugate, substrate, and stop solutions are ready to use and do not require any further dilutions.

4. Dilute the ABRAXIS® Wash Buffer (5X) Concentrate at a ratio of 1:5 with deionized or distilled water. If using the entire bottle (100 mL), add to 400 mL of deionized or distilled water and mix thoroughly.

F. Working Scheme

The microtiter plate consists of 12 strips of 8 wells, which can be used individually for the test. The standards must be run with each test. Never use the values of standards which have been determined in a test performed previously.

| | 1 | 2 | 3 | 4 | 5 | 6 | 7 | 8 | 9 | 10 | 11 | 12 |
|---|-------|-------|--------|--------|---|---|---|---|---|----|----|----|
| A | Std 0 | Std 1 | Std 5 | Sample | | | | | | | | |
| B | Std 0 | Std 1 | Std 5 | Sample | | | | | | | | |
| C | Std 1 | Std 5 | Sample | | | | | | | | | |
| D | Std 1 | Std 5 | Sample | | | | | | | | | |
| E | Std 0 | Cont. | etc. | | | | | | | | | |
| F | Std 0 | Cont. | etc. | | | | | | | | | |
| G | Std 0 | LRB | | | | | | | | | | |
| H | Std 0 | LRB | | | | | | | | | | |

Std 0 - Std 5: Standards

Contr.: Control (QCS)

LRB: Laboratory Reagent Blank

Samp1, Samp2, etc: Samples

G. Assay Procedure

1. **Add 50 µL of the standard solutions, control, LRB, or samples** into the wells of the test strips according to the working scheme given. Analysis in duplicate or triplicate is recommended.
2. **Add 50 µL of the antibody solution** to the individual wells successively using a multi-channel pipette or a stepping pipette. Cover the wells with parafilm or tape and mix the contents by moving the strip holder in a circular motion on the benchtop for 30 seconds. Be careful not to spill the contents. **Incubate the strips for 90 minutes at room temperature.**
3. Remove the covering, decant the contents of the wells into a sink, and blot the inverted plate on a stack of paper towels. **Wash the strips three times using the diluted wash buffer.** Please use at least a volume of 250 µL of 1X wash buffer for each well and each washing step. Blot the inverted plate after each wash step on a stack of paper towels. After the last wash/blot, check the wells for any remaining buffer in the wells, and if necessary, remove by additional blotting.
4. **Add 100 µL of the conjugate solution** to the individual wells successively using a multi-channel pipette or a stepping pipette. Cover the wells with parafilm or tape and mix the contents by moving the strip holder in a circular motion on the benchtop for 30 seconds. Be careful not to spill the contents. **Incubate the strips for 30 minutes at room temperature.**
5. Remove the covering, decant the contents of the wells into a sink, and blot the inverted plate on a stack of paper towels. **Wash the strips three times using the diluted wash buffer.** Please use at least a volume of 250 µL of 1X wash buffer for each well and each washing step. Blot the inverted plate after each wash step on a stack of paper towels. After the last wash/blot, check the wells for any remaining buffer in the wells, and if necessary, remove by additional blotting.
6. **Add 100 µL of substrate (color) solution** to the individual wells successively using a multi-channel pipette or a stepping pipette. Cover the wells with parafilm or tape and mix the contents by moving the strip holder in a circular motion on the benchtop for 30 seconds. Be careful not to spill the contents. **Incubate the strips for 20-30 minutes at room temperature.** Protect the strips from sunlight.
7. **Add 50 µL of stop solution** to the wells in the same sequence as for the substrate (color) solution using a multi-channel pipette or a stepping pipette.
8. Read the absorbance at 450 nm using a microplate ELISA photometer within 15 minutes after the addition of the stopping solution.

H. Evaluation

The evaluation of the ELISA can be performed using commercial ELISA evaluation programs such as 4-Parameter (preferred) or Logit/Log. For a manual evaluation, calculate the mean absorbance value for each of the standards. Calculate the %B/B₀ for each standard by dividing the mean absorbance value for each standard by the Zero Standard (Standard 0) mean absorbance. Construct a standard curve by plotting the %B/B₀ for each standard on the vertical linear (y) axis versus the corresponding Microcystins concentration on the horizontal logarithmic (x) axis on graph paper. %B/B₀ for the control (QCS), LRB, and samples will then yield levels in ppb of Microcystins by interpolation using the standard curve. Results can also be determined using a spreadsheet macro available from Gold Standard Diagnostics upon request.

The concentrations of the samples are determined using the standard curve run with each test. Samples showing a lower concentration of Microcystins than standard 1 (0.05 ppb) should be reported as containing < 0.05 ppb of Microcystins. Samples showing a higher concentration than standard 5 (5.0 ppb) must be diluted to obtain accurate results. The concentration of the positive control (QCS) provided should be 0.75 ± 0.185 ppb.

Semi-quantitative results can be derived by simple comparison of the sample absorbances to the absorbances of the standards. Samples with lower absorbances than a standard will have concentrations of Microcystins greater than that standard. Samples which have higher absorbances than a standard will have concentrations of Microcystins less than that standard.